

## The MELAG Euroklav® Series

The cost-effective »Class S« Autoclaves







"The reliable MELAG device saves considerable time and ensures the highest levels of safety for patients and practice team alike."





"We set the focus on integral health and the well-being of our patients. The MELAG System Solution perfectly supports this approach."

28CliniC,

"In MELAG, we have found an excellent partner providing the best-possible hygiene standards and optimal workflows."

Sparkle Dental Boutique, Abu Dhabi

"The MELAG System ensures a quick and safe decontamination workflow in our clinic."

> Prosthodontic Associates, Australia

## The Euroklav® Series

### Unique in price and performance.

The quick and reliable Euroklav® Series represents an even more affordable alternative for the cost-effective steam sterilization. Euroklav® autoclaves are thus the ideal solution for practices and clinics using instruments that do not require sterilization in a Class B device. The MELAG Class S autoclaves are equipped with a high-performance vacuum pump in order to ensure best sterilization and drying results of even wrapped instruments with simple hollow bodies.

Rely on our innovative solutions for instrument decontamination. On the following page, you will find the key advantages at a glance:






Trust:
Over 500,000 units sold
worldwide





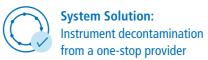



### The product highlights:

#### **Performance**

Quick sterilization and drying times enable efficient workflows.

### Safety

The integrated documentation interface provides maximum process reliability.

### **Flexibility**

Customize the chamber depths and the water supply to your individual needs.

#### **User-friendliness**

The MELAG 4-key concept and our digital tools make for an intuitive operating concept.



Find out more about the product highlights in our video: www.melag.com/en/multimedia



# Aligned to your workflow:

#### Euroklav® 23 VS+

Stand-alone steam sterilizer with a 45 cm chamber depth



#### Euroklav® 29 VS+

Stand-alone steam sterilizer with a 35 cm chamber depth





### Save even more time and money:

With an optional MELA*dem*® water treatment system, you benefit from time- and cost-saving operation through the automated water supply and disposal.

### The Euroklav® product highlights These advantages are unique.

Gain more time for your patients: All our research and development efforts are focused on providing an efficient and safe workflow in instrument decontamination. Benefit from the time- and cost-saving advantages provided by short sterilization times, the intuitive operating and the versatile documentation concept of our Class S autoclaves.



Short sterilization processes for your practice and clinic: Quick-Program (including pre-vacuum) sterilizes unwrapped instruments in only 15 minutes\*. Larger quantities of wrapped instruments can be reprocessed in only 25 minutes\* using the Universal-Program.

MELAG provides a wide range of insert racks to ensure an individual use of the sterilization chamber. The aligned accessories make for the optimal loading of the steam sterilizers with up to 4 kg of instruments.

\* without drying

MELAG steam sterilizers monitor all the relevant process parameters using a highly precise microprocessor control to guarantee highest safety standards. The program logs provide legal certainty and are saved via an serial interface.

Our aligned documentation media offer maximum flexibility: Customize the documentation process with MELAflash CF-Card Printer, MELAnet Box, MELAprint® 44 or the approval software MELAtrace® to your needs.

One system – two autoclaves: To meet the individual needs in sterilizing instruments, MELAG offers the Euroklav® in two sterilization chamber sizes: Euroklav® 29 VS+ with 17 litres and Euroklav® 23 VS+ with 23 litres.

The integrated water tanks and patented air-cooling system make the stand-alone steam sterilizer completely independent of installation requirements and easy to install. An additional MELAG water treatment unit for automated water supply reduces your workload and thus saves time and money.

The easy use of Euroklav® steam sterilizers allows you to concentrate on patient treatment: The proven MELAG **4-key concept** enhances your workflow and avoids operating

In addition to the intuitive operating concept, our **digital tools** take the instrument decontamination process to the next level. Discover these innovative multimedia solutions on the MELAG website.

### **Digital tools for** smart workflows

Benefit from the digital tools on our website.





Learn from more than 300 videos about safe instrument decontamination: www.melag.com/en/multimedia





## Download Center:






# The MELAG System Solution Your one-stop provider for instrument decontamination.

The Euroklav® steam sterilizers are an important element of an efficient instrument decontamination workflow. The perfect alignment of the MELAG System Solution helps you to save time and money whilst providing the best-possible protection for patients and staff.

Benefit from the MELAG System Solution with products from the leading manufacturer in infection control:





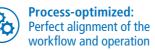

Enhance the quality of instrument decontamination: Combine the MELAG steam sterilizers with MELAtherm® washer-disinfectors for efficient cleaning and disinfection, MELAseal® sealing devices for safe instrument packaging and the user-friendly documentation and approval software MELAtrace® to implement an ideal workflow in practice and clinic.



### The instrument decontamination process:



Cleaning and disinfection with MELAtherm®



Packaging with MELAseal®



Sterilization with Euroklav®



Documentation and approval with MELAtrace®



Labelling with MELAprint® 60\*

### Perfection in the finest detail

### Configure your Euroklav® steam sterilizer.

Use the full potential of the Euroklav® Series: With the wide range of accessories available for the MELAG autoclaves, you can tailor the workflow to your requirements, thereby ensuring the best instrument decontamination workflow.





### **Loading system**

Take advantage of our loading system to ensure the optimal use of the sterilization chamber.



#### **Mount A Plus**

accommodates five trays for single

Art. no. 82620 for Euroklav® 29 VS+ Art. no. 82630 for Euroklav® 23 VS+



#### **Mount D**

can be loaded with two higher sterilization containers

Art. no. 46840 for Euroklav® 29 VS+ and Euroklav® 23 VS+



#### **Sterilization containers**

are perforated and equipped with paper filters for perfect sterilization and drying results. Discover our wide range of sterilization containers:

www.melag.com/en/products/ packaging/sterilization-containers



#### Foil holder

permits the vertical sterilization of wrapped instruments for optimal drying results

Art. no. 22410 for Euroklav® 29 VS+ Art. no. 22420 for Euroklav® 23 VS+



### **Documentation**

The Euroklav® steam sterilizers quarantee legally-compliant documentation and approval with modern output media.



#### MELA*trace*®

for the professional documentation, approval and traceability of the entire decontamination process

MELAtrace® without load administration: Art. no. 21128 MELAtrace® Pro with load administration: Art. no. 21138



#### MELAnet Box

transfers all logs to the computer and into the clinic network automatically

Art. no. 40296



#### MELAflash CF-Card Printer

saves all logs of the steam sterilizer simply and easily on a CF-Card Art. no. 01039



#### MELAprint® 44

prints the sterilization logs in paper form

Art. no. 01144



### Water treatment

MELAdem® water treatment units provide high-quality demineralized water to protect both your steam sterilizer and your instruments.



#### MELAdem® 40

is an ion exchanger for practices and clinics sterilizing up to three cycles per day

MELA*dem*® 40: Art. no 01049 Connection set: Art. no. 09031



#### MELAdem® 47

is a high-performance reverse osmosis unit for practices and clinics with larger sterilization requirements

MELA*dem*® 47: Art. no. 01047 Connection set: Art. no. 09032



#### MELA*dem*<sup>®</sup> 53 / 53 C

is a high-capacity water treatment unit for MELAtherm® and up to three further devices

Art. no. 01038



#### MELA*jet*®

facilitates instrument rinsing and flushing with demineralized water\* Art. no. 27300

## **Technical data**

### Figures and facts at a glance.

Programs		Euroklav® 23 VS+ Euroklav® 29 VS+	Drying time
Quick-Program S*	134°C, Holding time: 3,5 min.	15 min.	10 min.
Universal-Program	134°C, Holding time: 5,5 min.	25 min.	25 min.
Gentle-Program	121°C, Holding time: 20,5 min.	40 min.	25 min.

<sup>\*</sup>For unwrapped instruments only

	Euroklav® 23 VS+	Euroklav® 29 VS+	
Article number	10056	10057	
Device type	Stand-alone (additional water treatment unit can be connected)		
Chamber size	45 cm long x Ø 25 cm	35 cm long x Ø 25 cm	
Chamber volume*	23 litres	17 litres	
Load quantity	4 kg instruments 1 kg textiles	3 kg instruments 0,8 kg textiles	
Device dimensions (W x H x D)	$42.5 \times 49 \times 70 \text{ cm}$ Device fits on a 60 cm table top	$42,5 \times 49 \times 62$ cm Device fits on a 50 cm table top	
Weight	45 kg	42 kg	
Power supply	220-240 V / 50/60 Hz		
Power consumption	2.300 Watt	2.100 Watt	



PREMED PHARMA KFT.

CÍM 2040 Budaörs, Gyár u. 2.

TELEFON 06 23 889 700

FAX 06 23 889 710

E-MAIL info@premedpharma.hu

WEB www.premedpharma.hu





For further information: www.melag.com

MELAG Medizintechnik